
Care Connect Set-up

1. If you are not already connected to the Care Connect EMR application, you may download the desktop application via this link:


<https://services-acc.careconnect.be/client/6.2-gp/>

Download this zip file and extract all from it. Ensure that the extracted folder is called careconnect and is sitting on your Desktop.


2. Once you have logged in you will come to your home page.

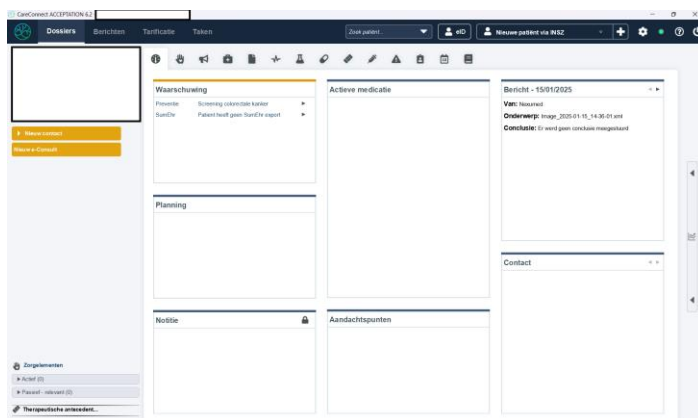3. Click the settings icon in the top right corner.

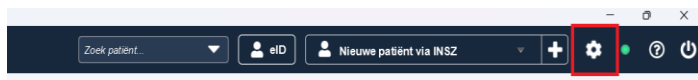4. You will see a drop down menu. Look for the titles Berichten and Interfaces.

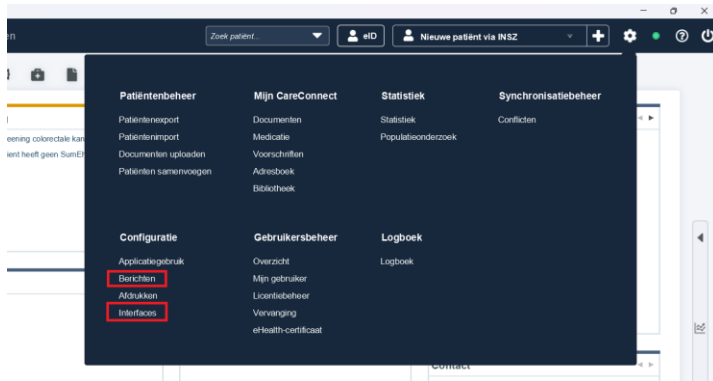5. First click Berichten and you will see this displayed on the page but Mexi interface will be a lighter colour. Click it as seen in the photo at which time the font will become bold.

6. 1. click the first(1)input field and copy and paste this:

C:\Nexumed\inFromDevice\output-to-emr

2. click the second (2) input field and copy and paste the same thing

C:\Nexumed\inFromDevice\output-to-emr

3. The third (3) red square will then be a bold colour and click that to save this address's

Now Berichten will be properly configured. Next, we must configure the Interfaces.

7. Click the setting wheel/icon again to get the menu dropdown. This time click the Interface from the dropdown. You come across this on the screen.

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Click the checkbox Gebruik unidirectionele toestelkoppeling as since in the red box above.

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Cardioflex
ECG F12
 SDS104
 SpiroConnect
 Winspiro Pro

9. This pop up will then be displayed on the screen

Toestelkoppeling:

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Software argumenten

1. Click the first (1) input field and copy and paste:


C: path to your Desktop => C:\Nexumed\emrReciever\dummy.exe

2. Click the second (2) input field and copy and paste:

C:\Nexumed\nexumedIn\

3. Last click the third red box (3) and click OK

10. The pop up will disappear and you will see this and click the red box to save the interface configuration.

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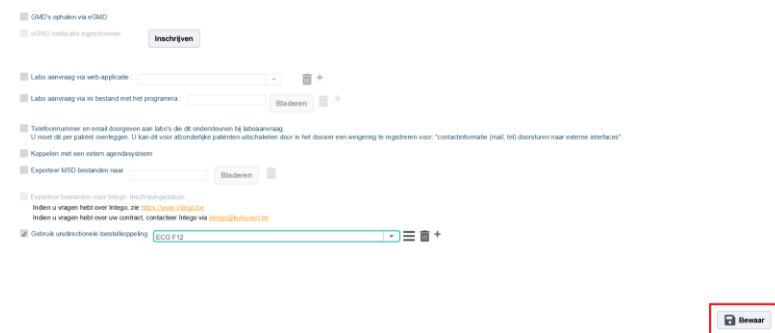☒ Druk de ondersteunde beschikking

Opslaan

This will take you to the end of the configuration of your Care Connect desktop application with Nexcore software by Nexumed.

Enjoy

The Nexumed Team

www.nexumed.eu